

ARTWORK FILE

PROJECT NAME

CAVEX

DOCUMENT NAME

IMPRESAFE_IFU HOLDER OUTSIDE

VERSION

2.1


DATE

130622

PRINT TYPE

TBD

COLOURS

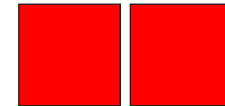

PMS 289 PMS SILVER

MATT FINISH PAPER

BARCODE

-

FOR INFORMATION ONLY
NOT FOR PRODUCTION

NAAF ONT WERP

© NAAF ONTWERP 2022

CAVEX IMPRESAFE

-
Effective spectrum:
Bactericidal | Yeasticidal | Virucidal*

Complies with EN 13624, EN 13727, EN 14476,
EN 14561, EN 14562

Composition in 100 g: 30.0 g Benzalkonium chloride,
auxiliaries and water

Exposure time: 3 % solution in 3 minutes



United Disinfectant Manufacturers AG
Allmendstrasse 21 8320 Fehraltorf
Switzerland +41 (0)55 511 03 33
info@udm.swiss www.udm.swiss



United Disinfectant Manufacturers AG
Dr. Grass-Strasse 12, 9490 Vaduz,
Principality of Liechtenstein
hello@udm.li www.udm.li



Manufactured exclusively for



Cavex Holland BV
Fustweg 5 2031 CJ Haarlem
The Netherlands +31 23 530 77 00
dental@cavex.nl cavex.nl



1 liter

* Limited virucidal
(incl. HIV/HCV/HBV)



see bottle



2460

CAVEX

HA0501.07



cavex.nl



2460



CAVEX IMPRESAFE

-
UFI: KYW-28MC-100F-TSC7

DANGER

H302 Harmful if swallowed.

H314 Causes severe skin burns and
eye damage.

H410 Very toxic to aquatic life with
long lasting effects.

P280 Wear protective gloves/protective
clothing/eye protection/face
protection.

P301+P312 IF SWALLOWED: Call a POISON
CENTER or doctor/physician if
you feel unwell.

P301+P330 IF SWALLOWED: rinse mouth.

+P331 Do NOT induce vomiting.

P303+P361 IF ON SKIN (or hair):

+P353 Take off immediately all
contaminated clothing. Rinse skin
with water/shower.

P305+P351 IF IN EYES: Rinse cautiously
with water for several minutes.

+P338 Remove contact lenses, if
present and easy to do.
Continue rinsing.

P501 Dispose of contents/container in
accordance with local and
national regulations.

ARTWORK FILE

PROJECT NAME

CAVEX

DOCUMENT NAME

IMPRESAFE_IFU HOLDER INSIDE

VERSION

2.1


DATE

130622

PRINT TYPE

TBD

COLOURS

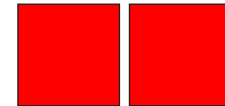

PMS 289

PMS SILVER

MATT FINISH PAPER

BARCODE

-

FOR INFORMATION ONLY
NOT FOR PRODUCTION

NAAF ONT WERP

© NAAF ONTWERP 2022

CAVEX IMPRESAFE

-

Instructions

1. Squeeze bottle to dispense 30 ml
2. Pour 970 ml of water (always first) in reservoir
3. Pour disinfectant into reservoir
4. Place impression in basket
5. Lower the basket with the impression into the fluid
6. Leave impression for 3 minutes in the fluid for complete disinfection
7. Take impression out of fluid and rinse thoroughly under running water of at least drinking quality for 10 seconds
8. Shake off excess water and put impression in a plastic bag. Seal the bag airtight and mark the bag as "disinfected"



FOR PROFESSIONAL USE ONLY.
Not for use on invasive medical devices. Do not mix with other products. Make sure that all visible soiling has been removed before disinfection.



Store in cool, dry and well-ventilated facility. Keep out from direct sunlight, heat, and open flame. Storage temperature: <30°

A list of compatible and incompatible materials is available on our website. If you are uncertain as to the material composition of the medical device, contact the manufacturer before use.