

CAVEX IMPRESAFE

Effective spectrum:

Bactericidal | Yeasticidal |
Virucidal against enveloped viruses
Complies with EN 13624, EN 13727, EN 14561, EN 14562, EN 17111

Exposure time: 3% solution in 3 minutes

Composition in 100 g: 30.0 g benzalkonium chloride, auxiliaries and water

DANGER

- H302 Harmful if swallowed.
- H314 Causes severe skin burns and eye damage.
- H410 Very toxic to aquatic life with long lasting effects.
- P280 Wear protective gloves/protective clothing/eye protection/face protection.
- P301+P312 IF SWALLOWED: Call a POISON CENTER or doctor/physician if you feel unwell.
- P301+P330 +P331 IF SWALLOWED: Rinse mouth. Do NOT induce vomiting.
- P303+P361 +P353 IF ON SKIN (or hair): Take off immediately all contaminated clothing. Rinse skin with water/shower.
- P305+P351 +P338 IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing.
- P501 Dispose of contents/container in accordance with local and national regulations.

Manufactured exclusively for

 Cavex Holland BV • Fustweg 5 • 2031 CJ Haarlem
The Netherlands
+31 23 530 77 00 • dental@cavexnl • cavex.nl

CAVEX IMPRESAFE

- EN DANGER
- AR خطر
- BG ОПАЧО
- CS NEBEZPEČÍ
- DA FARE
- DE GEFÄHR
- EL ΚΙΝΔΥΝΟΣ
- ES PELIGRO
- ET ETTEVAATUST
- FI VAARA
- FR DANGER
- HR OPASNOST
- HU VESZÉLY
- IT PERICOLO
- KO 위험
- LT PAVOJINGA

- LV BRIESMAS
- NL GEVAAR
- NO FARE
- PL NIEBEZPIECZEŃSTWO
- PT PERIGO
- RO PERICOL
- RU ОПАЧО
- SK NEBEZPEČENSTVO
- SL NEVARNO
- SR OPASNOST
- SV FARA
- TH ภัยอันตราย
- TR TEHLIKE
- UK НЕБЕЗПЕКИ
- ZH-CN 危险
- ZH-TW 危險



UFI: KYYW-28MC-100F-TSC7



GTIN 10000000000001
YYYY-MM-DD
LOT 00-123456




United Disinfectant Manufacturers AG
Dr. Grass-Strasse 12 • 8490 Vaduz
Principality of Liechtenstein
Tel. +423 237 15 03 • info@udm.li • www.udm.li




An affiliate of United Disinfectant Manufacturers AG
Allmendstrasse 21 • 8320 Fehraltorf • Switzerland
Tel. +41 (0)55 511 03 33 • info@udm.swiss
www.udm.swiss

HA050i.11


Instructions

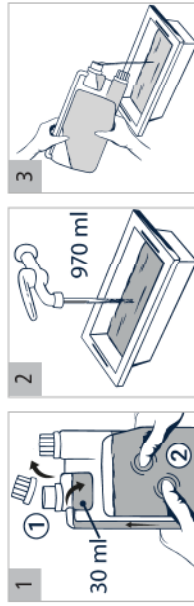


1 Squeeze bottle to dispense 30 ml.

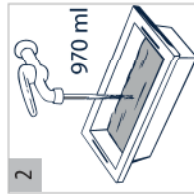


2 Pour 970 ml of water (always first) in reservoir.

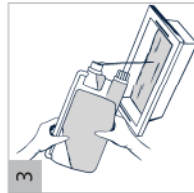


3 Pour disinfectant into reservoir.

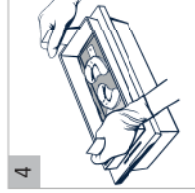

4 Place impression in basket. Lower the basket with the impression into the fluid.

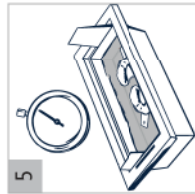

5 Leave impression in the fluid for 3 minutes for complete disinfection.

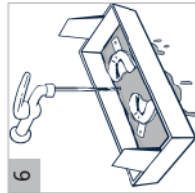

6 Take impression out of fluid and rinse thoroughly under running water of at least drinking-water quality for 10 seconds.



7 Shake off excess water and put impression in a plastic bag. Seal the bag airtight and mark the bag as "disinfected".

Refer to the instructions for use for more information.



FOR PROFESSIONAL USE ONLY.

Do not mix with other products. Ensure that all visible soiling has been removed before disinfection. Failure to adhere to these limitations may reduce the efficacy of the product, thereby increasing the risk of infection.



Store in a cool, dry and well-ventilated place. Keep away from direct sunlight, heat, and open flame. Storage temperature: < 30 °C. Dispose of the product after expiry date.

DOSING TABLE

3%	30 ml	+	970 ml		3 min
----	-------	---	--------	--	-------

A list of compatible and incompatible materials is available on cavex.nl. If you are unsure about the material composition of the medical device to be disinfected, contact the manufacturer before use.

